



Achieve Korean Glass Skin **with GENOSYS approach**

**This manual should be used only for training of GENOSYS users.

**Revised in November 2025

Any part of this user manual must not be modified, reproduced, or copied without prior permission of DTS MG Co., Ltd.
(Hereafter referred to as DTS MG)

Content

1. What is Korean Glass Skin and how to achieve?
 2. How Professional Microneedling Helps Achieve Glass Skin?
 3. How Professional No-Needle Microneedling Treatment Can Improve Skin?
 4. Importance of Homecare
 5. Importance of Safe Formulation for Achieving Korean Glass Skin
 6. Exfoliate without Over-Cleansing
 7. Moisturize During the Day and Night
 8. Protect, Nourish and Regenerate
 9. Why Step-by-Step SPF and Coverage Products Are Beneficial?
-

1. What is Korean Glass Skin and how to achieve?

“Korean Glass Skin” is a popular beauty ideal that originated from South Korea’s skincare culture. It describes skin that looks clear, luminous, smooth, and almost translucent — like glass.

This look is not achieved with makeup, but through deep hydration, consistent exfoliation, and strong skin barrier health. The goal is to make the skin so healthy and plump that it naturally reflects light, giving that dewy, radiant finish.

- Available treatments which help achieve Korean Glass Skin

The Korean glass skin look is achieved through consistent skincare and advanced aesthetic treatments that focus on hydration, barrier repair, and collagen renewal.

- **Water Peel - HydraFacial or Aqua Peel**
This treatment deeply cleanses the pores using water-based suction and exfoliation. It infuses the skin with hydrating serums and antioxidants, leaving it clear, smooth, and dewy immediately after the procedure.
 - **Skin Boosters (Hydration Injections)**
Examples include Rejuran Healer, Profhilo, Juvelook, NCTF, or hyaluronic acid boosters. These are injected microdroplets of polynucleotides or hyaluronic acid that improve moisture retention, elasticity, and glow. Over time, they stimulate collagen and enhance skin density and translucency.
 - **Laser Toning or Picotoning**
Lasers such as Q-switched Nd:YAG, PicoSure, or Picocare target pigmentation, dullness, and uneven tone. They refine skin texture, tighten pores, and promote collagen renewal, helping to create a clearer and more reflective complexion.
 - **Microneedling with Serums**
Microneedling treatments performed with microneedling devices or microneedle rollers (GENOSYS) are often combined with hyaluronic acid, peptides, or growth factor serums. This promotes skin regeneration, minimizes pores, and improves absorption of active ingredients, resulting in plump and radiant skin.
 - **Fractional Laser (Non-ablative or Sub-ablative)**
Examples include Fraxel Dual, ResurFX, and Laveen Thulium Laser. These treatments improve texture, fine lines, and uneven tone while boosting dermal collagen to create smoother, more light-reflective skin.
 - **Cryotherapy or Oxygen Infusion**
Oxygen therapy delivers pressurized oxygen combined with nutrients or hyaluronic acid mist, while cryotherapy tightens pores and soothes inflammation. Both give the skin a fresh and glazed appearance.
 - **LED or Low-Level Light Therapy**
Red or near-infrared LED light enhances cellular regeneration and circulation. It boosts overall radiance and calms redness after other procedures.
 - **Custom Korean Combination Treatments**
Many Korean clinics combine different methods, such as laser toning with Rejuran Healer and Aqua Peel, or microneedling with exosome infusion and LED therapy.
-

These approaches aim to create comprehensive improvements in skin clarity, texture, and glow.

- At-home Maintenance

Consistent at-home care is essential to maintain the glass skin effect. This includes multi-layer hydration using essences, serums, ampoules, masks and creams, along with daily sun protection.

2. How Professional Microneedling helps achieve Korean Glass Skin?

Microneedling is a minimally invasive clinical procedure that uses fine, sterile needles to create controlled micro-injuries in the skin. These tiny punctures trigger the skin's natural healing response and stimulate collagen and elastin production. Over time, this leads to smoother, firmer, and more even skin texture.

Reasons microneedling supports the glass skin effect:

- ✓ Refines skin texture
 - Reduces the look of pores, acne scars, and fine lines, giving the surface a smoother, light-reflecting quality.
- ✓ Boosts collagen and elasticity
 - Encourages fibroblast activity to produce new collagen and elastin, helping skin appear firmer and more supple.
- ✓ Enhances serum absorption
 - The microchannels created during treatment allow hydrating and repairing ingredients like hyaluronic acid, peptides, or growth factors to penetrate deeper.
- ✓ Improves skin tone
 - Helps fade uneven pigmentation and post-acne marks, resulting in a clearer, more luminous complexion.
- ✓ Strengthens the skin barrier
 - When followed with soothing, repairing serums, microneedling supports long-term barrier health, not just short-term glow.

GENOSYS has a combination of the best microneedle tools, post care tools and specialized cosmeceuticals which are formulated to be used with microneedling procedures.




How the active ingredient of GENOSYS Power Solutions – sh-polypeptide 7 can bring Korean glass skin effect with microneedling?

Sh-polypeptide-7, a bioengineered peptide that mimics human growth hormone, promotes skin regeneration and collagen synthesis. When used together with microneedling, it can significantly enhance the rejuvenation process and help achieve the smooth, luminous look associated with Korean glass skin.

Microneedling creates controlled microchannels that stimulate the skin's natural repair response, and these channels allow sh-polypeptide-7 to penetrate deeply into the epidermis and dermis. This deeper delivery increases its bioavailability and amplifies its effects on fibroblast activity, collagen production, and overall skin renewal.

The combination supports improved dermal density and firmness, leading to plumper, more elastic skin. It also refines texture, reducing the appearance of pores and fine lines, and helps the surface reflect light evenly, giving the skin a glass-like clarity.

Overall, the synergy between sh-polypeptide-7 and microneedling produces both immediate luminosity and long-term rejuvenation, creating firmer, smoother, and more translucent skin that embodies the Korean glass skin aesthetic.

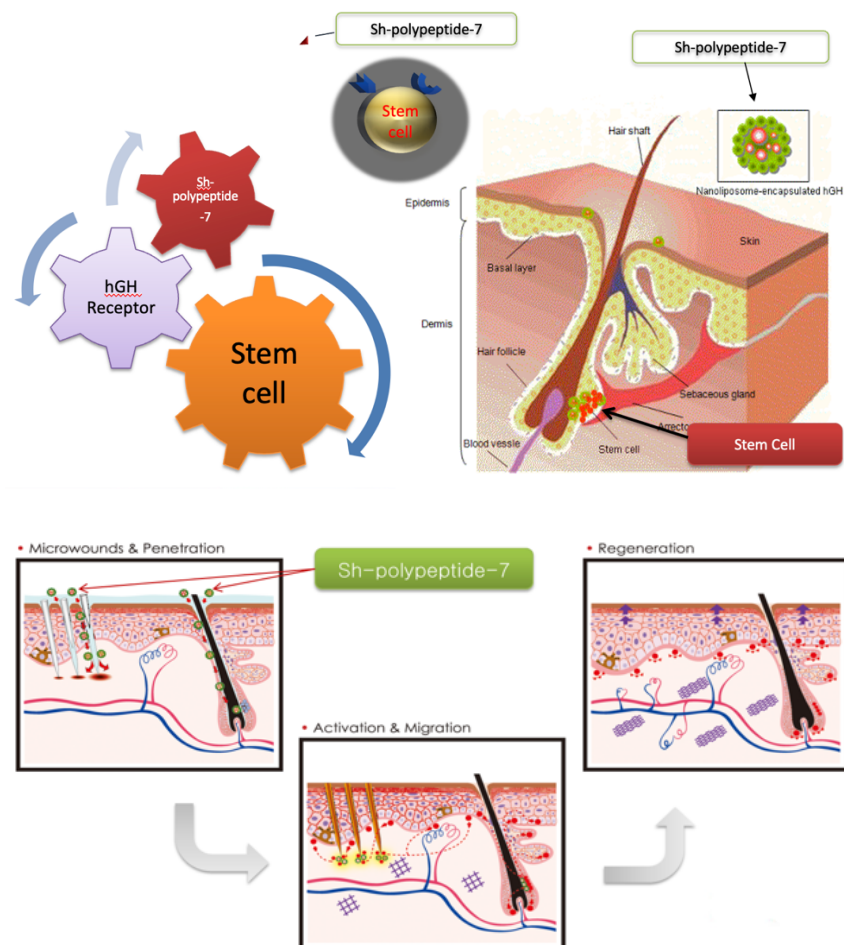

<How sh-polypeptide 7 works in the skin>

Professional GENOSYS Microneedling Treatment Step

Step 1. SKIN DEFENDER LIP & EYE MAKEUP REMOVER

Step 2. SNOW O₂

Step 3. SNOW BOOSTER / INTENSIVE PROBLEM CONTROL TONER

Step 4. EZ CO₂ MASK – Pre-needling carboxy treatment

Option. Microneedling with SRS – Fractional peeling

Step 5. Microneedling with POWER SOLUTIONS – AWS. HES. SWS. PCS / CTS.CVS – Cocktail method

Step 6. PEPTIDE GEL MASK with GENO-LED IR 2 - Post needling treatment

Step 7. SOOTHING REPAIR POST CREAM / CREAMS

Step 8. MULTI SUN CREAM / ULTRA SHIELD SUN CREAM

Step 9. INTENSIVE BLEMISH BALM CREAM

Post care - MICROBIOME ENERGY INFUSING MIST(DURING THE DAY) – Prevention of TEWL

Different Professional Treatment Options

Basic Microneedling Treatment	Roller/Needle Pen K + 1 Solution (AWS/SWS/PCS) Peptide Gel Mask / Soothing Bomb Sea Algae Mask Soothing Repair Postcream Multi Sun Cream / Ultra Shield Sun Cream
Cocktail Microneedling Treatment	Roller/Needle Pen K + 2 Solutions (AWS/SWS/PCS/HES + CTS/CVS) Peptide Gel Mask / Soothing Bomb Sea Algae Mask Soothing Repair Postcream Multi Sun Cream / Ultra Shield Sun Cream
Carboxy Cocktail Microneedling Treatment	EZ CO ₂ Mask Roller/Needle Pen K + 2 Solutions (AWS/SWS/PCS/HES + CTS/CVS) Peptide Gel Mask / Soothing Bomb Sea Algae Mask Soothing Repair Postcream Multi Sun Cream / Ultra Shield Sun Cream

**Carboxy Acid Cocktail
Microneedling
Treatment**

EZ CO2 Mask
SRS Peeling
Roller/Needle Pen K + 2 Solutions (AWS/SWS/PCS/HES + CTS/ CVS)
Peptide Gel Mask / Soothing Bomb Sea Algae Mask
Soothing Repair Postcream
Multi Sun Cream / Ultra Shield Sun Cream

<GENOSYS Professional Treatment product range>
Pre Needling Products

Products for Needling

Products for Bio Needling

Post Needling Products

3. How professional no-needle microneedling treatment can improve skin?

A professional spicule treatment (no-needle microneedling treatment) supports skin's renewal process by accelerating cell turnover, enhancing product absorption, and improving skin texture without harsh exfoliation or downtime.

Spicules are microscopic needle-like structures derived from freshwater sponges. When applied to the skin, they create thousands of microchannels that lightly stimulate the dermis — similar to microneedling, but using a topical method instead of a device.

These spicules remain in the skin temporarily (usually 24–72 hours) and work from within to trigger regeneration and collagen production.

GENOSYS offers a specialized professional treatment, enriched with the uniquely formulated BIO-MESO™ PDRN(60000ppm), panthenol(10000ppm), ceramides and an anti-aging complex to regenerate skin and strengthen the skin barrier.

- How GENOSYS BIO-MESO™ PDRN treatment creates the Glass Skin effect

- ✓ Stimulates skin renewal
Spicules gently exfoliate the outer layer while activating cellular turnover in the deeper layers. This removes dull surface cells and promotes the growth of fresh, new skin — resulting in a smoother, clearer surface that reflects light more evenly.
- ✓ Boosts collagen and elasticity
The controlled micro-stimulation from spicules encourages collagen synthesis. Over time, this strengthens the dermal structure, improving firmness and reducing fine lines for a tighter, “lifted” glass-like finish.
- ✓ Enhances absorption of active ingredients
The microchannels created by spicules allow serums and ampoules — such as hyaluronic acid, peptides, or niacinamide — to penetrate more deeply, maximizing hydration and nourishment. Deep hydration is key to achieving that luminous, light-reflective glow.
- ✓ Refines pores and texture
By promoting renewal and tightening the skin, spicule therapy minimizes the appearance of enlarged pores and uneven texture, giving the skin a porcelain-smooth appearance.
- ✓ Improves tone and brightness
Regular spicule treatments help reduce pigmentation and uneven tone by accelerating the natural shedding of melanin-rich cells, revealing a clearer and more radiant complexion.

***** Tip : Treating the skin with BIO-MESO™ PDRN Homecare Ampoule 5000 can enhance the professional treatment effects.



<Professional Bio-Meso PDRN>



<Homecare Bio-Meso PDRN>

4. Importance of Homecare

Professional treatments provide powerful results, but homecare maintains and enhances them. A consistent at-home skincare routine, inspired by Korean layering techniques, ensures that the skin remains hydrated, balanced, and glowing between sessions.

Essential homecare steps include:

- ✓ Gentle and effective cleansing – to remove sunscreen and impurities without stripping the skin
- ✓ Hydrating toner or essence – to replenish moisture and prepare the skin for active ingredients.
- ✓ Gentle exfoliation – to gently exfoliate the skin for boosting cell turnover
- ✓ Serum or ampoule – containing hydrating and repairing ingredients such as hyaluronic acid, growth factors, niacinamide, or peptides.
- ✓ Moisturizer / nourishing – to lock in hydration with a lightweight but nourishing cream.
- ✓ Sunscreen – to prevent pigmentation and protect the skin's clarity.

Consistent homecare helps preserve the collagen and hydration developed during microneedling sessions, ensuring long-term radiance and skin strength.

<GENOSYS Homecare range>

Cleansing & Peeling				
				
DMR	SOC	SBT	PCT	EPG
Treatment Program - Serum & Cream				
				
Anti-Wrinkle	Skin Depigmenting	Hydrating	Calming & Nourishing	Anti acne & Combination skin

Special Eyecare



Sunscreen & Coverage



Mist & Mask



GENOSYS Morning Routine

- Step 1. SNOW O₂
- Step 2. SNOW BOOSTER / INTENSIVE PROBLEM CONTROL TONER
- Step 3. SERUMS – AWS/MVS/PCS/AFS/MHS
- Step 4. EYECCELL EYE COUNTOUR SERUM
- Step 5. EYECCELL EYE CONTOUR CREAM
- Step 6. CREAMS – MFC / MVC / PCC / HSC / MHC / EOC
- Step 7. NDCell ANTI-WRINKLE CREAM on Neck and Décolleté
- Step 8. MULTI SUN CREAM / ULTRA SHIELD SUN CREAM on Face & Neck
- Step 9. INTENSIVE BLEMISH BALM CREAM
- Step 10. SKIN CARING BLEMISH BALM CUSHION
- Step 11. MICROBIOME ENERGY INFUSING MIST(DURING THE DAY)

GENOSYS Evening Routine

- Step 1. SKIN DEFENDER LIP & EYE MAKEUP REMOVER
- Step 2. SNOW O₂
- Step 3. EPI TURNOVER BOOSTING PEELING GEL
- Step 4. SNOW BOOSTER / INTENSIVE PROBLEM CONTROL TONER
- Step 5. MASK – PEPTIDE GEL MASK / SOOTHING BOMB SEA ALGAE MASK / BIO-FERMENT AGE DEFYING POWDER MASK
- Step 6. SERUMS – AWS/MVS/PCS/AFS/MHS
- Step 7. EYECCELL EYE COUNTOUR SERUM
- Step 8. EYECCELL EYE CONTOUR CREAM
- Step 9. CREAMS – MFC / MVC / PCC / HSC / MHC / EOC
- Step 10. NDCell ANTI-WRINKLE CREAM on Neck and Décolleté
- Step 11. SKIN RESCUE OVERNIGHT CREAM MASK

5. Importance of safe formulation for achieving Korean Glass Skin

The foundation of the Korean glass skin look is healthy, balanced, and resilient skin. True radiance does not come only from treatments—it comes from skin that is strong, calm, and free from chronic irritation.

That is why safe formulations—free from parabens, ethanol, artificial fragrance, artificial pigment, and sulfates—are essential. GENOSYS products keep having safe formulations whenever possible. They help maintain the skin's natural barrier, prevent inflammation, and ensure long-term clarity and glow.

1) Protecting the skin barrier

The skin barrier is one of the most important factors in achieving a smooth, light-reflective surface. Harsh preservatives or surfactants can weaken this barrier, leading to dryness, redness, and sensitivity.

By using clean formulations, the skin's natural lipid layer stays intact, keeping moisture inside and irritants out. This is the key to maintaining the dewy, hydrated "glass skin" appearance.

2) Preventing irritation and inflammation

Ingredients such as ethanol (alcohol), artificial fragrance, or pigment can trigger irritation or micro-inflammation, especially in sensitive or post-procedure skin. Chronic irritation results in uneven tone, redness, and loss of translucency.

Gentle, fragrance-free products allow the skin to remain calm, even-toned, and more receptive to active ingredients that promote glow and firmness.

3) Supporting skin health after professional treatments


After microneedling or spicule treatments, the skin becomes more permeable and sensitive. Using products without potential irritants ensures faster healing and safer absorption of beneficial nutrients.

Clean formulations reduce the risk of post-treatment reactions and help maintain the smooth, luminous results achieved through professional care.

4) Ensuring long-term radiance

The glass skin effect is not a temporary shine; it reflects balanced skin health. Long-term use of safe, gentle skincare helps stabilize the skin's microbiome, prevents sensitivity buildup, and supports consistent hydration and elasticity.

This results in skin that naturally glows without relying on heavy makeup or filters.

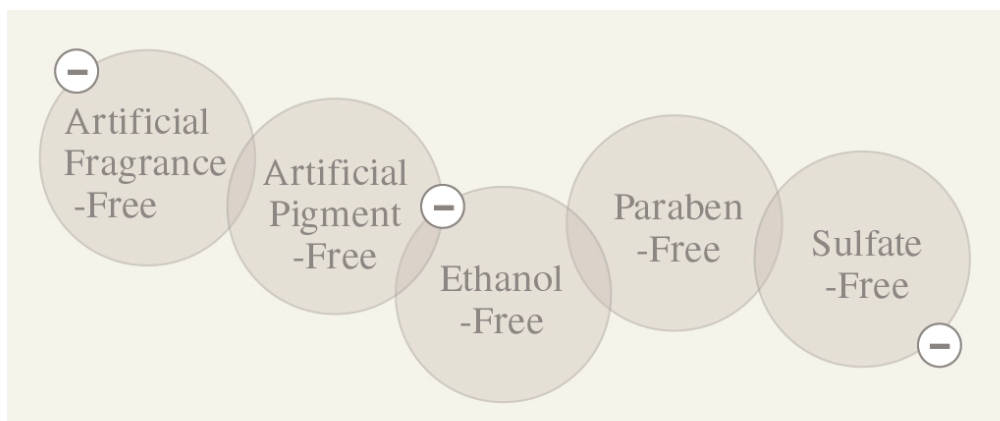

<5 No additions concept of GENOSYS formulation>


6. Exfoliate well without over-cleansing

Exfoliating without over-cleansing helps maintain the skin's natural barrier and hydration. Over-cleansing removes essential lipids and disrupts the microbiome, which can lead to dryness, sensitivity, and dullness. Gentle exfoliation clears away dead cells while preserving the skin's protective oils. Using a mild cleanser once or twice a day is enough to keep the skin clean. Exfoliating 1 to 2 times a week with enzyme-based exfoliants keeps the skin smooth and radiant. Harsh scrubbing or cleansing until the skin feels tight should be avoided. After exfoliation, applying hydrating toners or essences helps restore moisture. This balance between cleansing and hydration allows the skin to stay clear, healthy, and naturally luminous.

GENOSYS has an excellent exfoliating product. GENOSYS EPI TURNOVER BOOSTING GEL is an enzyme peeling gel that removes dead skin cells without irritation. It contains moringa, so-called 'Miracle Tree' which purifies and nourishes skin and the plant complex from desert moisturizes and soothes skin.

Benefits :

- Peeling gel without rough particles → No discomfort when scrubbing
- Moist gel type without leaving dryness
- Removing dead skin cells.
- Binding and removing waste materials

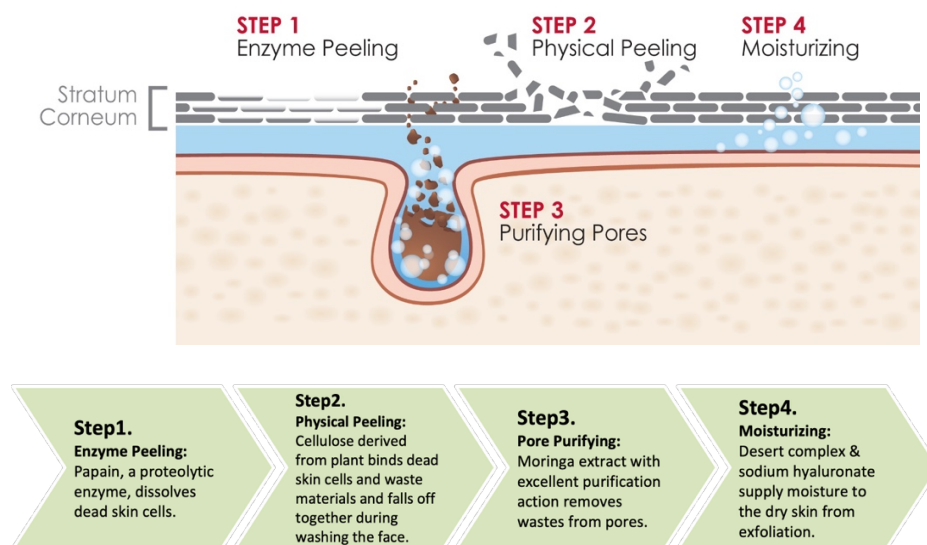

<Mechanism of EPI TURNOVER BOOSTING GEL>



< GENOSYS EPI TURNOVER BOOSTING GEL>

7. Moisturize during the day and night

Moisturizing the skin both during the day and at night is essential for maintaining the smooth, hydrated, and radiant quality that defines glass skin. The skin constantly loses moisture through evaporation and exposure to environmental stress, so consistent hydration keeps the barrier strong and light-reflective.

During the day, moisturizer protects the skin from dehydration caused by sunlight, pollution, air conditioning, and temperature changes. It helps create a smooth surface that reflects light evenly, giving the skin a fresh and luminous appearance. A day moisturizer also acts as a base for sunscreen and makeup, preventing dullness and flaking throughout the day.

At night, moisturizing supports the skin's natural repair and renewal process. While we sleep, the skin becomes more permeable and active in regenerating cells, so applying a rich or barrier-repairing moisturizer helps restore lost lipids and strengthen the structure of the outer layer. This nighttime hydration prevents dryness, improves elasticity, and ensures that the skin remains plump and glowing by morning.

Keeping the skin consistently moisturized both day and night allows it to stay balanced, smooth, and radiant, creating the translucent, light-reflective appearance that characterizes Korean glass skin.

i. During the day - MICROBIOME ENERGY INFUSING MIST

Using a microbiome mist during the day helps maintain the balance of healthy bacteria on the skin, which is essential for keeping the barrier strong and hydrated. A balanced skin microbiome supports natural defense mechanisms, reduces inflammation, and helps the skin retain moisture more effectively. Throughout the day, the skin is exposed to environmental stressors such as pollution, UV radiation, and temperature changes, all of which can disturb its microbial balance and cause dryness or dullness. Spraying a microbiome mist replenishes hydration without disrupting makeup or sunscreen, while also restoring the skin's protective ecosystem. Regular use keeps the surface smooth, supple, and resilient, which enhances light reflection and contributes to the clear, dewy appearance.

GENOSYS MICROBIOME ENERGY INFUSING MIST revitalizes skin and enhances skin's natural strength with unique blend of pre/probiotics that corrects the balance of skin microbiome and hyaluronic acid complex that replenishes the skin's moisture barrier.

Effects :

- Powerful moisturizing effect
- Strengthens and protects skin barrier

- Soothes skin and restores skin condition
- Gives skin natural glow and radiance



<GENOSYS MICROBIOME ENERGY INFUSING MIST>

ii. During the night - SKIN RESCUE OVERNIGHT CREAM MASK

SKIN RESCUE OVERNIGHT CREAM MASK containing pink ceramide, growth factors, pumpkin fruit extract, centella asiatica, and phytosphingosine works synergistically to repair, hydrate, and rejuvenate the skin overnight—key elements for achieving glass skin.

Pink ceramide helps restore the lipid barrier, locking in moisture and improving elasticity so the skin stays smooth and supple. Growth factors stimulate cell renewal and collagen synthesis, supporting firmness and a more refined texture. Pumpkin fruit extract provides natural enzymes and vitamins that gently exfoliate dead cells while nourishing the skin with antioxidants, promoting a brighter, more even tone. Centella asiatica soothes inflammation, strengthens the skin's resilience, and supports healing, reducing sensitivity and redness. Phytosphingosine, a natural component of the skin's barrier, helps regulate microbial balance and further reinforces protection against environmental stress.

Together, these ingredients maintain deep hydration, strengthen the barrier, and enhance skin clarity overnight. The result is a fresh, plump, and luminous complexion by morning, reflecting light evenly and achieving the translucent, healthy glow.



<GENOSYS SKIN RESCUE OVERNIGHT CREAM MASK>

8. Protect, nourish and regenerate

Products that regenerate, protect, and nourish the skin are vital for achieving glass skin because they build lasting skin health and radiance. Regenerative ingredients promote cell renewal and collagen production, creating a smooth, light-reflective surface. Protective formulas strengthen the barrier and shield against environmental damage, preventing dryness and irritation. Nourishing products supply essential lipids and nutrients that maintain hydration and elasticity. Together, these products keep the skin balanced, resilient, and luminous.

GENOSYS offers various protective and regenerative products.

1) SOOTHING REPAIR POSTCREAM

GENOSYS SOOTHING REPAIR POSTCREAM is a regenerating cream that promotes skin recovery after professional treatment. It helps irritated skin rapidly recover from redness, erythema and edema after the dermatological operations or any other irritations and promotes healthy rejuvenation.

Features

- A finishing cream of invasive medical aesthetic treatments
- Sh-Polypeptide-7 rapidly penetrates deep into skin and activates potential inactive stem cells
- Excellent effect in anti-inflammation and anti-oxidation
- Helps fast restoration of damaged skin tissue
- Helps irritated skin rapidly recover from redness, erythema and edema after treatment



<GENOSYS SOOTHING REPAIR POSTCREAM>

2) EGF REPAIR OXYMASK CREAM

GENOSYS EGF REPAIR OXYMASK CREAM is a unique bubbling mask cream which makes dull & stressed skin healthy with skin regenerating ingredients such as EGF and madecassoside. Oxygen is capable of improving cellular metabolism, accelerating the healing process, reducing skin irritations, and producing an anti-inflammatory effect.

Benefits :

- A cream which regenerates skin damaged by various cause (S.O.S cream)
- Supplying oxygen to skin
- Skin regenerating
- Soothing
- Moisturizing

SOOTHING REPAIR POSTCREAM	EGF REPAIR OXYMASK CREAM
Intensive repair for a week right after the professional treatment	Daily regenerating cream for the skin damaged by various causes (S.O.S cream)



< GENOSYS EGF REPAIR OXYMASK CREAM>

3) SKIN BARRIER PROTECTING CREAM

GENOSYS SKIN BARRIER PROTECTING CREAM encourages healthy and soft skin by retaining moisture and protecting skin barrier with enriched ceramide, botanical extracts and amino acid complex.

Benefits :

- Promotes water retention
- Supplies nutrients
- Anti-inflammatory effect
- Normalizes the intercellular lipids



<GENOSYS SKIN BARRIER PROTECTING CREAM>



<Clinical test results of GENOSYS SKIN BARRIER PROTECTING CREAM>

9. Why Step-by-Step SPF and coverage products are beneficial?

The Korean skincare philosophy focuses on layering lightweight products rather than relying on a single heavy formula. This approach ensures better protection, natural coverage, and a flawless glow — essential for achieving the “glass skin” look.

GENOSYS has a 3-step SPF system — Sunscreen(Multi Sun Cream/ Ultrashield Sun Cream) , BB Cream(Intensive Blemish Balm Cream), and BB Cushion(Skin Caring BB Cushion) — works together to protect, perfect, and maintain the skin’s radiant finish throughout the day.

1st step : Sunscreen (SPF 40–50) – The Skin Protector

GENOSYS MULTI SUN CREAM (SPF 40 PA ++) & **ULTRA SHIELD SUN CREAM (SPF 50+ PA ++++)**

Sunscreen is the first and most crucial layer. A broad-spectrum SPF shields the skin from harmful UVA and UVB rays that cause premature aging, pigmentation, and collagen breakdown.

Benefits:

- Creates the base layer of defense against photoaging
- Prevents dark spots and dullness
- Keeps skin smooth and even-toned, which is essential for maintaining the glass skin clarity
- GENOSYS Suncreams have light texture and use chemical+mineral filters
- When applied as the first layer, it ensures full protection before makeup or tinted products, GENOSYS Suncreams do not create a white cast effect

GENOSYS MULTI SUN CREAM (SPF 40 PA ++)

Sun protection product which gives a smooth, non-greasy application sensation with excellent adhesiveness to the skin

- Effective sun protection
- Paraben-Free, Alcohol-Free, Pigment Free
- Soothing effect
- Makes skin cells healthy after UV damages
- Can be used without irritation after medical treatment. (microneedling, laser, peeling treatment)



<GENOSYS MULTI SUN CREAM>

GENOSYS ULTRA SHIELD SUN CREAM (SPF 50+ PA ++++)

Sun protection product which strongly defends skin against UV rays and improves sunburn

- Broad spectrum sun protection
- Sunburn care effect
- Reef-safe: not containing harmful ingredients that affect coral bleaching → Suitable for swimming and marine sports

- Can be used without irritation after medical treatment (microneedling, laser, peeling treatment)



<GENOSYS ULTRA SHIELD SUN CREAM>

2nd step : BB Cream – The Skin Perfector

GENOSYS INTENSIVE BLEMISH BALM CREAM (SPF30 PA++)

BB Cream bridges skincare and makeup. GENOSYS INTENSIVE BLEMISH BALM CREAM has a skin caring formula and a specially blended colour has been added to this formula. It provides light coverage, tone correction, and added SPF, while still delivering hydration and skin-nourishing benefits.

Benefits:

- Evens out skin tone and blurs minor imperfections
- Adds a second layer of UV protection
- Contains beneficial ingredients like adenosine, Phytalex SC, allantoin and oregano leaf extract
- Keeps the skin breathable and natural-looking, unlike heavy foundations that can clog pores
- Self adapting to own skin colour
- Hydration – minimizing moisture loss
- Skin Protection from harmful environment



<GENOSYS INTENSIVE BLEMISH BALM CREAM>

3rd step : BB Cushion – The Skin Refresher

GENOSYS SKIN CARING BLEMISH BALM CUSHION (SPF 50+ PA ++++)

The BB cushion provides convenient, on-the-go coverage and sun protection touch-ups. It contains SPF and hydrating ingredients, making it ideal for reapplying protection throughout the day without disturbing the base layers.

Benefits:

- Refreshes makeup while reapplying SPF, which is a key K-beauty secret for consistent sun protection
- Delivers a dewy, luminous finish that maintains the glass skin glow
- Lightweight and easy to use for quick touch-ups
- Adds a soft-focus filter effect for a poreless, smooth look
- GENOSYS SKIN CARING BLEMISH BALM CUSHION is a BB cushion that can be used after professional treatment. More than 60% of the product is composed of moisture essence, which enables a natural and healthy glow. Various peptide complex helps calm the irritated skin.
- When used on top of the GENOSYS INTENSIVE BLEMISH BALM CREAM, a long lasting protection and coverup effect is maximized





<GENOSYS SKIN CARING BLEMISH BALM CUSHION>



<GENOSYS multi-layer SPF and cover-up system>